

*Replaceable Cartridge Abator*

## Scrub your emissions 95% free of ethylene oxide gas!

We are pleased to introduce the newest accessory for your Andersen sterilizer.



The Replaceable Cartridge Abator is designed to remove more than 95% of the ethylene oxide emitted from 200 AN79 Anprolene ampoules. Installation is simple, as the abator attaches easily to most new **AN74i, AN74ix, AN74j, AN2000** or **EOGas 4** model sterilizers.

The emissions from ethylene oxide are regulated in many parts of the world. If these regulations apply to your area, the Replaceable Cartridge Abator may be the right solution for you.

---

### Installation:

The vent hose from the sterilizer connects to the AN5100 Abator which is vented to the outside of the building for extra safety.

Cabinet dimensions: 12"W x 26"H x 17"L


### **AN5100 Abator**

**Order #: 115v - AN5100.01**

**230v - AN5100.11**

### **Replacement Cartridge:**

**Order #: AN5100.00**

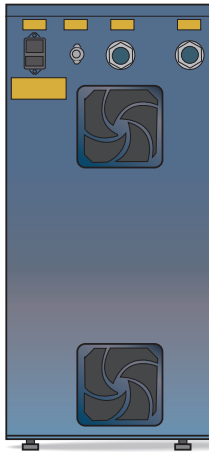


### Installation is Quick & Easy

The AN5100 Replaceable Cartridge Abator is controlled and monitored through a single wire connection to the sterilizer. Accurate tracking of gas usage will automatically alert the operator when it approaches time to change the cartridge.



### Disposable Abator Cartridge

The AN5100.00 Replacement Cartridge is available for shipping by normal means. After 200 cycles the old cartridge can be disposed of with ordinary waste if local regulations permit.

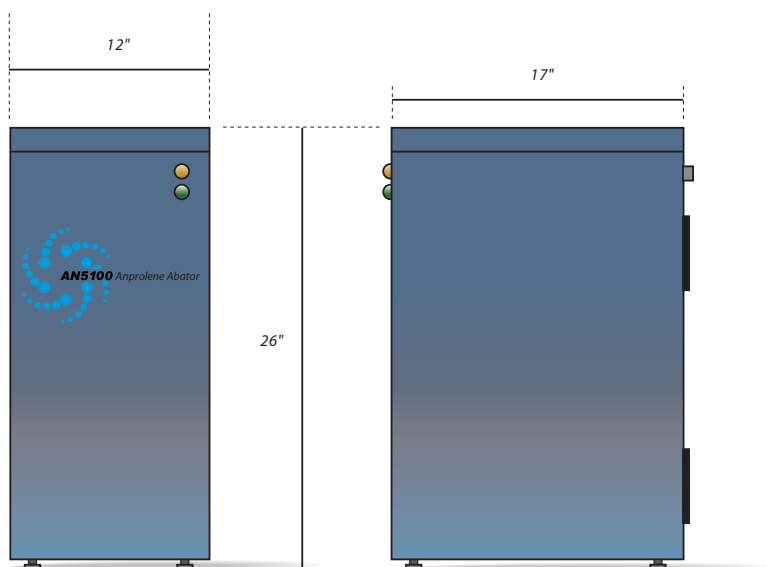

#### Specifications:

115V, 2 amps

Indoor use, temperature range 68 °F - 91 °F

Operating weight 67 lbs.

AN5100.00 Replacement Cartridge

Manufactured by: Andersen Sterilizers, Inc.

*AN74i/ix units starting with lot number 243588/243583, 9/14/2004 are abator ready and include the BNC connector for the AN5100 Replacement Cartridge Abator.*

*All AN74j and AN2000 and EOGas 4 Sterilizers are abator ready.*

**Anprolene**<sup>®</sup>  
gas sterilization

**EOGas**<sup>®</sup> 4

PN667\_071207