



It's soft, gentle,
and shorter



Try the world's softest Nasogastric tube. Difficult to insert? Ask for the Andersen tube with stylet.

- The new, shorter length AN11.5 facilitates care of patients in incubators.
- The standard AN11 Andersen Pediatric Tube 4 ft. length has been shortened to eliminate tangling and kinking.
- Double lumen (vented), marked at 20, 30 and 40 cm from the most proximal, aspirating port, radiopaque.
- All Andersen Tubes come equipped with an anti-reflux filter.

New Shorter Length

This new Andersen Tube was designed for neonatal patients and infants where, in an incubator, there is little room for extra tubing. The new, shorter length is 24 inches (60 cm). This helps to eliminate tangling in small spaces and promote better management of nasogastric drainage.


It is marked at 20, 30, and 40 cm from the most proximal aspirating port to permit more accurate placement in infants.

Sizes

The original Andersen Tube is 48 inches (120 cm) and is available in French sizes 10 and 16. The Andersen Tubes include a Stylet in Adult, 16 French (AN10S) and Pediatric, 10 French (AN11S).

Double Lumen Construction

All are vented (double lumen), radiopaque and come equipped with an anti-reflux filter. Best known for the soft PVC tubing, the Andersen Tube's vent lumen is exterior to the aspirating lumen, a more physiological shape, yielding substantially more aspirating lumen than another brand of the same exterior size. The aspirating ports act to screen out particles which are too big to pass through the aspirating lumen.

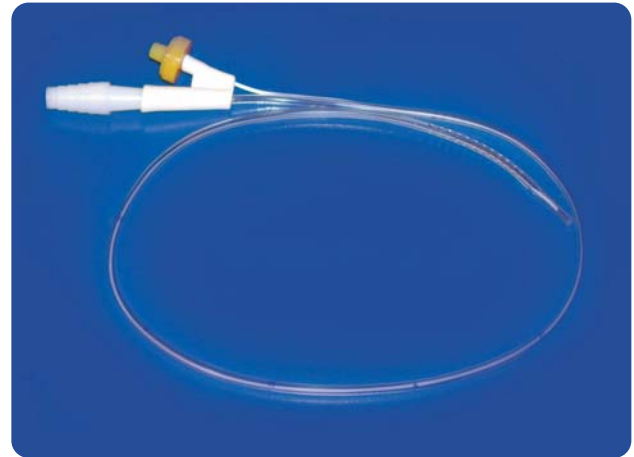

Need a STYLET?

Andersen has the only nasogastric tube with a Stylet to improve intubation in the unconscious patient. It is stiff when you pass it. Remove the stylet and it is soft when the patient awakens. (See sizes)

Anti-Reflux Filter

The air vent on all Andersen Tubes comes equipped with a unique non-wettable filter that acts to prevent tube reflux.

Reliability

Visual assurance of the Andersen Tube function is simple. The need for routine irrigation, and routine checks for tube patency or pump function are eliminated. "If it's bubbling, it's working!"

REORDER No.	SIZE	DESCRIPTION
AN10	16Fr.	The most efficient size for most adult patients.
AN10S	16Fr.w/Stylet	
AN11	10Fr.	Pediatric, with marks at 20, 30, 40, and 50 cm.
AN11S	10Fr.w/Stylet	
AN11.5	10Fr. Extra Short	Pediatric, with marks at 20, 30, and 40.

www.anpro.com

customerservice@anpro.com

 **ANDERSEN
PRODUCTS**
the future of gas sterilization

www.anpro.com
Andersen Products, Inc.
Health Science Park • 3202 Caroline Drive
Haw River, NC 27258-9564 USA
800-523-1276 • fax 336-376-8153

www.andcal.com
H. W. Andersen Products of California, Inc.
Health Science Park • 3151 Caroline Drive
Haw River, NC 27258-9575 USA
800-524-3455 • fax 336-376-3088

Distributed in Europe by:
H. W. Andersen Products, Ltd.
Davy Road • Clacton-on-Sea
Essex CO15 4XA UK